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#### This unit must be used for laboratory applications only and for internal use only.

The manufacturer declines all responsibility for any use of the unit that does not comply with these instructions.

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#### This unit has been designed and manufactured in compliance with the following standards:

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UL 61010-1

CAN/CSA-C22.2 No.61010-1

VELP reserves the right to modify the characteristics of its products to constantly improve their quality.

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## 1. Introduction

With VELP TEMPSoft™ the user can constantly control and manage the temperature, through the wireless interface. Many aspects can be monitored, including the set point control (max. and min. temperature alarm thresholds), the monitoring of the internal temperature trend, the possibility to set working ramps with different temperatures and times and store test data in conformity with GLP.

The communication between the software and the instrument is managed by the Wireless DataBox™.

## 2. Assembling and installing

Upon receipt and after having removed the packaging, please check the integrity of the instrument. The box includes:

- Wireless DataBox™
- USB cable
- TEMPSoft™

The Wireless DataBox™ comes already prepared and doesn't need any action for the assembling it. Using the USB cable is possible connect it to the PC. On the USB-Pen Drive is possible to find the driver if required by Windows after the first connection to the PC.

## 3. TEMPSoft™

### 3.1 Installation

General rules for installing the TEMPSoft™ Software:

In order to install the TEMPSoft™ Software correctly, administrator rights for the PC must be available. The PC on which the TEMPSoft™ Software will be installed must not be overloaded with other software. This is to prevent unstable operation of the TEMPSoft™ Software.

To install the TEMPSoft™ Software:

1. Place the USB-Pen Drive in the USB-Port.
2. To start the installation run **setup.exe** from the USB-Pen drive.
3. The installation guide is displayed. Click **Next** and follow the instructions. The TEMPSoft™ icon will be created on the screen automatically.
4. Connect the Wireless DataBox™ to the PC.
5. After a few seconds the PC acknowledges the presence of a new hardware and will automatically start the installation procedure of the drivers. When requested, specify the path of the **Driver** folder on the USB-pen drive and proceed with installation following the instructions that appear on the screen.
6. Re-boot the computer to put the settings provided by the installation procedure into effect.

### 3.2 Turning on

Before turning on the instrument, all the steps described in the previous paragraphs must be carried out:

1. Turn on the PC and the screen.
2. Check that the TEMPSoft™ has been installed.
3. Connect the Wireless DataBox™ to the PC using the USB cable.
4. Launch the TEMPSoft™ (shortest way > double-click on the appropriate icon on the screen) (Fig. 1).
5. The "User" level password is "velp". Enter the password and press OK (Fig. 2).
6. The main window will appear (Fig. 3).

Fig. 1

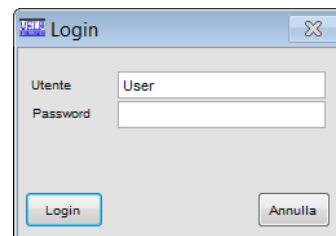
Fig. 2

With the level User, is not possible delete any analysis from database, or any temperature methods. In order to do this, is necessary logon with the "Superuser" level. The "Superuser" level password is "velp01".

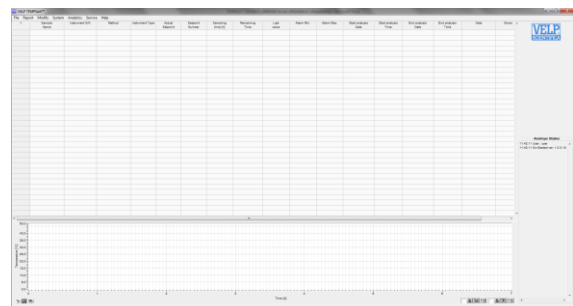
Fig. 3

## 4. TEMPSoft™ Quick Start Guide

This page shows an abstract of the TEMPSoft™.

The purpose is to help the user identify the various possibilities offered by the instrument and how to reach them easily.

<b>Menu</b>	<b>Submenu</b>	<b>Description</b>
File	New	Create a new analysis data base and "Laboratory information" file.
	Open	Load an already existing database.
	Save as	Save an analysis data base (including "Laboratory information" file) assigning a file name. If the file name already exists, it can be overwritten.
	Laboratory information	Enter the details of the laboratory concerned. This information will be included in the test report.
	Exit	Exit from the program.
Report	Create report	Create a standard report with the analyses selected.
	Export data	Export data of the rows selected to different format file.
Modify	Filter	Create database filters.
	Standard column	Reset the order of the database columns to the default setting.
	Configure	Configure the sequence of the database columns, show/hide columns.
	Remove	Remove one or more rows from the database.
	Select all	Select all the rows of the database containing a value.
System	Languages	Change the language of the program.
	RF Frequency	Permits to select the RF frequency (EU or US).
	ID Table	Permits to register a instrument S/N
Analytics	Methods	Display, modify, create or delete different temperature method.
Service	Connect Receiver	Manual search for Wireless DataBox™.
	Disconnect Receiver	Disconnect the Wireless DataBox™ from the PC.
	Acoustic signals	Select the operations that activate an acoustic signal.
	Login	Modify the software entry level.
	Diagnostic	Perform a diagnosis of the status of the instrument.
	Control panel	For VELP or authorized technical staff only.
Help	Software activation	Under construction.
	About	Contact VELP Scientifica directly by e-mail.
	Manual	PDF of TEMPSoft™ presentation.

## 5. DataBox™ and TEMPSoft™ Connection

In order to connect the Wireless DataBox™ to the PC, is necessary to use the USB cable, supplied with the unit.

1. After having connected it to the PC, select **Connect Receiver** from the **Service** menu (Fig. 4).
2. Using the "Connect Receiver" menu with the Automatic Search option, it is possible to find and connect automatically the Wireless DataBox™ (Fig. 5).
3. When the unit is connected to the TEMPSoft™, the software shows the S/N and a green led on the left bottom (Fig. 6).

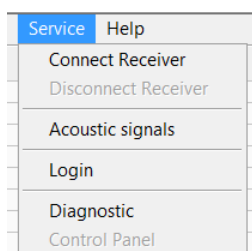
Fig. 4

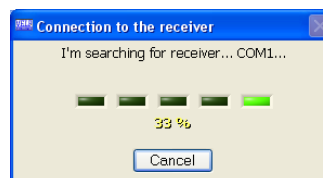
Fig. 5

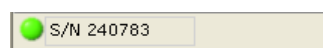
Fig. 6

## 6. TEMPSoft™ Workflow

### 6.1 Preparing for analysis

About the sample preparation, please see the analytical methods section of operating manual of main instrument.

#### 6.1.1 Filling the database

Several databases can be created and stored in the PC. Each database can store max 1000 analyses. Select **File/New** to create a new database or select **File/open** to load an already existing database.

Each database is connected to a “Laboratory information” file (Fig. 7). Before starting the analysis, some of the columns must be filled in. The following database information is mandatory and **must be entered** before starting the analysis:

<b>Sample Name</b>	Name of sample (alphanumeric field).
<b>ID</b>	ID of instrument that want to use with TEMPSoft™. Register S/N of instrument inside the ID Table.
<b>Method</b>	Name of the temperature test.
<b>Sampling Time</b>	Sampling time for the analysis.

Once the analysis has been carried out the content of the database columns cannot be modified. All the other database columns are filled automatically by the instrument at the end of analysis.

### 6.1.2 ID Instrument

The ID of the instrument that is necessary for connect to TEMPSoft™, must be set inside the ID Table, from the menu System. This value identifies uniquely an instrument.

### 6.1.3 Creating and choosing the Method

It is possible create different temperature’s method by using the Menu Analytics/Methods. In this menu is necessary indicate a name for the method, the type of instrument, the number of ramps and for any ramp, give the temperature and the time (Fig. 8).

Fig. 8

### 6.1.4 Choosing the “Sampling Time”

Using the **Sampling Time** it is possible to select different period of sampling, from 5 seconds to 24 hours (check software for all steps available). Fast sampling times require less testing in progress. For example 5 and 10 second permits only one test per time.

### 6.2 Performing analysis

In order to start the analysis, select the row and select START with the right button of mouse. Before performing an analysis, all the steps described above must be carried out (Fig. 9).

	Sample Name	ID	Method	Instrument	Sampling Time
1	Test Temperature	1	BOD 5	FOC120	10 sec

Fig. 9

The graph at the bottom of the display always refers to the analysis selected and not necessary to the analysis in progress.

### 6.2.1 Stopping analysis

In order to stop the analysis, select the row and with the right button of mouse, select STOP.

Once the analysis has been stopped, the icon  will be showed on the left part of analysis.

### 6.2.2 Alarms

TEMPSoft™ shows on screen a short explanation of the problem occurred in case of error.

### 6.3 Data evaluation

In order to evaluate the data obtained, the software offers the following possibilities:

- Configuring the data base
- Filtering the data base
- Exporting data

The indicator in the first column of the data base gives the status of each analysis in the series:

Indicator	Meaning
	Analysis ended correctly.
	Analysis terminated manually with a reset by user.
	Analysis with some notes.

### 6.3.1 Configuring the data base

The database consists of a set of data entered by the operator as well as data calculated by the software after analysis. Select **Modify** to modify the configuration of the database (Fig. 10).

<b>Filter</b>	Filter data.
<b>Standard column</b>	View columns in standard configuration.
<b>Configure</b>	Change the order of columns or conceal.
<b>Select all</b>	Selecting all the rows in the database.

Fig. 10

### 6.3.2 Filtering the database

The data shown can be filtered using the filtering function. Select **Modify/filters**. The following window appears (Fig. 11).

Select filters as follows:

1. Choose the fields to be filtered.
2. Select the operation to be performed on the field selected before and the value of reference.
3. If necessary, select AND/OR and choose the second option of filtering. (Fig. 12).
4. Press "Save New" in order to record the selection.
5. Select the new filter created or one already available.
6. To activate the filter selected Active Filter.

**NOTE:** to deactivate filter/s, select the filter/s as described above and press "Filter OFF".

The filter status ("Filter ON" or "Filter OFF") is shown at the bottom of the main window.

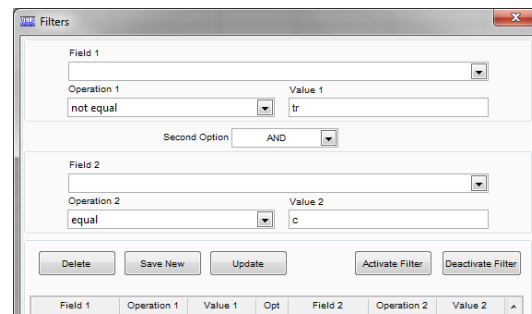
Fig. 11

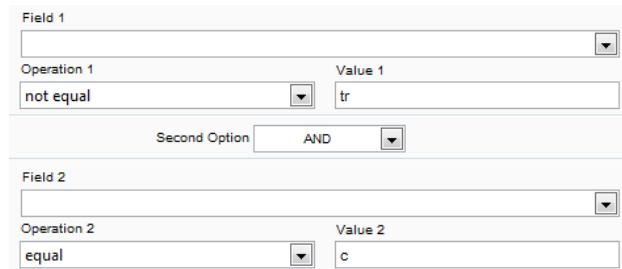
Fig. 12

### 6.3.3 Exporting data

Information can be exported from the data base in Excel, Text or Lims format by selecting **Report/Export Data EXCEL**, **Export Data TEXT** or **Export Data LIMS** (Fig. 13).

<b>Export Data EXCEL</b>	a .xls file is created; graphs not exported.
<b>Export Data TEXT</b>	a .txt file is created; graphs not exported.
<b>Export Data LIMS</b> (Fig. 17)	a .csv file format is created which can be easily read by the majority of laboratory management systems (LIMS). Select the value/unit to be exported (see below). Select the configuration of the exported data from the options available. It is possible to attribute the comma, semicolon or tab as delimiting character and the comma or decimal point as the decimal indicator depending on the type of the Lims. Press OK to confirm.

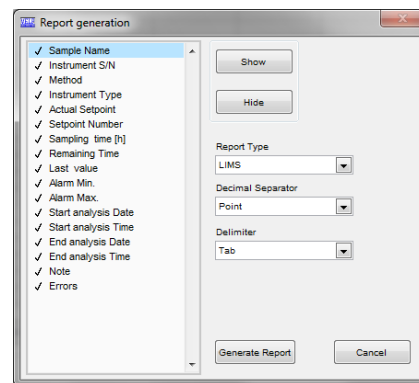
Fig. 13

## 7. Maintenance

The Wireless DataBox™ and TEMPSoft™ are totally maintenance-free.

## 8. Technical Data

Dimensions (WxHxD)	125x25x70 mm (4.9x1.0x2.7 in)	
Weight	85 g (0.04 lb)	
Power supply	5 V – USB Cable	
Power	0,225 W	
Overvoltage category	II	
Pollution degree	2	
Max. altitude	2000 m	
Max. humidity	85%	
Admitted room temperature	Storing -25 ÷ +65° C	Working +5 ÷ +40°C

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## 9. Parti di ricambio

10005112 Cavo USB A/Micro B 1.5m

We, the manufacturer VELP Scientifica, under our responsibility declare that the product is manufactured in conformity with the following standards:

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EN 61010-2-051/2003

EN61326-1/2013

and satisfies the essential requirements of the following directives:

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2012/19/UE (RAEE)

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## Thank you for having chosen VELP!

Established in 1983, VELP is today one of the world's leading manufacturer of analytical instruments and laboratory equipment that has made an impact on the world-wide market with Italian products renowned for innovation, design and premium connectivity. VELP works according to **ISO 9001**, **ISO14001** and **OHSAS 18001** Quality System Certification.

Our instruments are manufactured in Italy according to the IEC 1010-1 and CE regulation.

Our product lines:

### Analytical instruments

Elemental Analyzers  
Digestion Units  
Distillation Units  
Solvent Extractors  
Fiber Analyzers  
Dietary Fiber Analyzers  
Oxidation Stability Reactor  
Consumables

### Laboratory Equipment

Magnetic Stirrers  
Heating Magnetic Stirrers  
Heating Plates  
Overhead stirrers  
Vortex mixers  
Homogenizers  
COD Thermoreactors  
BOD and Respirometers  
Cooled Incubators  
Flocculators  
Overhead Shakers  
Turbidimeter  
Radiation Detector  
Open Circulating Baths  
Pumps

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VELP Scientifica Srl